





## **QuanTB Installation on a PC**

- 1. Visit the QuanTB website at <a href="http://www.msh.org/resources/quantb">http://www.msh.org/resources/quantb</a> to request a copy of QuanTB. You will receive an e-mail with a link to download the software. Use the link in the e-mail to download and save the QuanTB software to your computer.
- 2. Double-click on the QuanTB.zip file to unzip it. You will get a folder called QuanTB Download.
- 3. In the QuanTB Download folder, double-click on the **QuanTBforPC.exe** file to run it. The **Language Selection** screen will appear:



4. Select the language you want and click **OK** to continue. The following screen appears:











5. Click **Next** to continue then read the license agreement that appears. You must accept the terms of the agreement to continue the installation, as shown in the screen below.

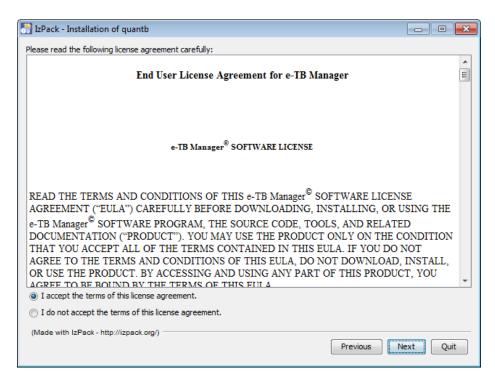

6. Enter an installation path for the software, as shown on the screen below and then click **Next**. We recommend accepting the default path that the program suggests.

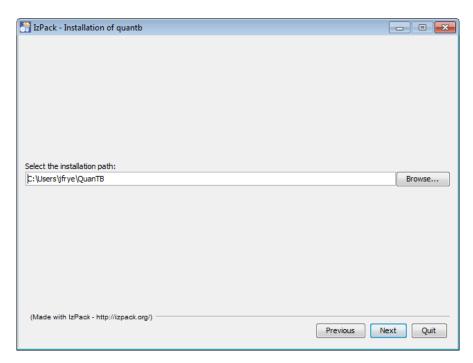








7. The program will confirm the installation directory it will create. Click **OK** to continue.

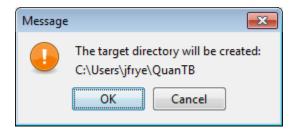


1. Select the packs you want to install. We recommend installing all of them. The greyed options are required, so you cannot unselect them. QuanTB comes with pre-populated lists of medicines and regimens standardized by the Global Drug Facility (GDF) and the World Health Organization (WHO) in the dictionary folder. QuanTB also has some "fictitious" forecasting for training/exercise purposes in the documents folder. If you have already installed QuanTB on your computer with customized lists of medicines and regimens or saved forecastings, your information will be overwritten if a new version is installed and the boxes are checked during the installation.

If you are reinstalling the software after customizing the medicines, regimens, or forecasting, we recommend that you uncheck the documents and dictionary boxes to preserve your existing files. Please read the Upgrading or Reinstalling QuanTB section of the user guide for more information. If this is the first time you are installing QuanTB, we recommend keeping the boxes checked. Click **Next** to continue.








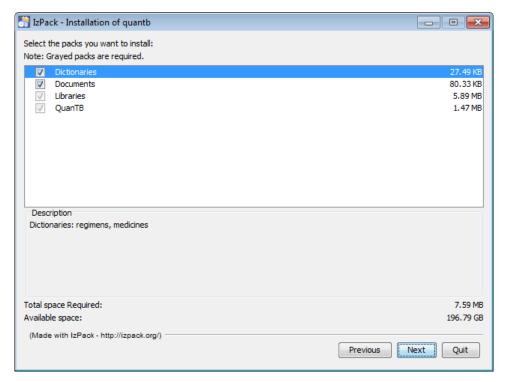

2. If you already have QuanTB installed in the same location, the system will warn you that you will overwrite files. Select No if you want to keep previously modified drug lists and/or treatment regimens.

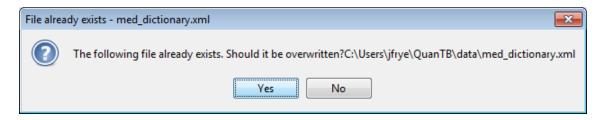


3. The software will display its installation progress, as shown on the screen below. Click **Next** to continue.







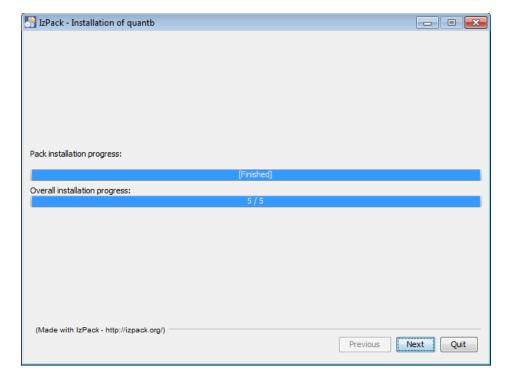









5. Choose which shortcuts you want the program to create for you, then click **Next**:

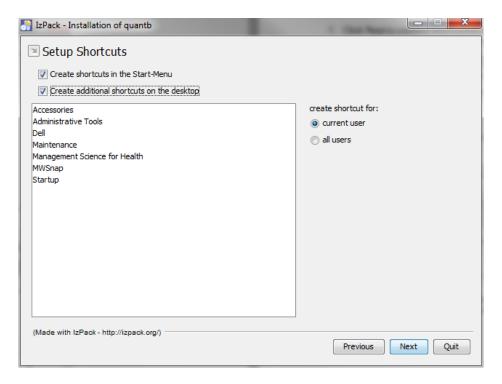

6. A screen confirming a successful installation will appear next.

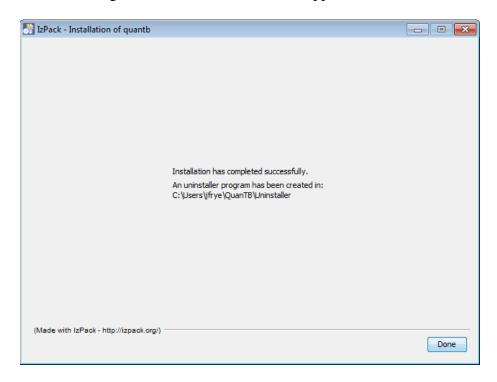







Depending on your version of Windows, you may get a Program Compatibility Assistant message after the installation (shown below). You should click the "This program installed correctly" option.

